

Traceable®

Traceable® Vaccine-Trac™ Data Logging Thermometer Bundle for Vaccines for Children (VFC) Program

Focus on high-quality care while ensuring VFC compliance

- Meets all CDC recommendations for digital data loggers in compliance with Vaccines for Children (VFC) program*
- Works with the TraceableLIVE® web platform to upload, manage and store reports
- Includes NIST-traceable certificate from our A2LA accredited (1750.01) calibration laboratory

*For more details, see the [Vaccine Storage and Handling Toolkit](https://www.cdc.gov/vaccine-storage-and-handling-toolkit/), CDC.gov (www.cdc.gov) is your online source for credible health information and is the official website of the Centers for Disease Control and Prevention (CDC).

thomassci.com



833.544.SHIP (7447)


Thomas
Scientific

Traceable® Vaccine-Trac™ Data Logging Thermometer Bundle for Vaccines for Children (VFC) program

This unique data logging bundle conforms to all CDC recommendations for vaccine storage and handling. (For more details, see the [Vaccine Storage and Handling Toolkit](#)). Bundle includes Vaccine-Trac data logger with bottle probe, bottle stand for optimal probe placement, USB drive for data transfer, and helpful signage to make compliance easier. The data logger comes with an NIST-traceable certificate from our A2LA accredited (1750.01) calibration laboratory. The logger and bottle probe are attached via a 10-ft (3-m) micro-thin cable so your refrigerator door can securely close without any issues.

Easily monitor temperatures overnight or for any time period—minimum and maximum temperatures are stored until the memory is cleared. Set high and low alarms in 1° increments; alarm sounds when temperature rises above or falls below set points. The alarm tracking feature stores data for up to 10 different alarm events detailing when the alarm state occurred, as well as when the unit returned within range. Alarm state indicator includes visual LEDs, audible alerts, and flashing LCD segments.

Transfer recorded data in a CSV format to any computer via a USB flash drive—no software is required. Simply clear memory after downloading data to USB drive. Data logger continues to record new data while downloading old data. Other status indicators include memory full, USB data transfer, active alarm state, and low battery.

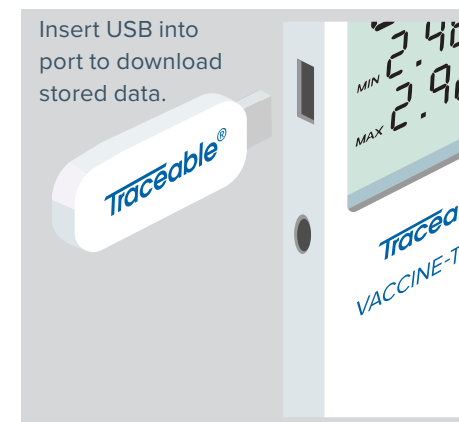

Ordering Information

Range	Resolution	Accuracy	Sampling rate	Memory capacity	Power	Item number
-58 to 158°F (-50 to 70°C)	0.01°	±0.25°C	5 seconds	525,600 points	Two AAA batteries and AC adapter (all included)	22A00G251

The CDC and the Vaccines for Children Program (VFC) are not affiliated with Traceable, and they do not evaluate or endorse the products sold on this website. CDC.gov (www.cdc.gov) is your online source for credible health information and is the official website of the Centers for Disease Control and Prevention (CDC).



ThomasSci.com 833.544.SHIP (7447)
CustomerService@ThomasSci.com

Connect With Us:

