

CDC16/APC6 Polyclonal Antibody

Catalog number: 14307-1-AP

Size: 20 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:


CDC16; ANAPC6, APC6, CDC16,

CDC16 homolog, CDC16Hs,


Cyclosome subunit 6

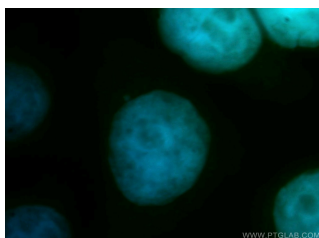

Immunofluorescent analysis of HeLa cells, using CDC16 antibody 14307-1-AP at 1:50 dilution and FITC-labeled donkey anti-rabbit IgG (green). Blue pseudocolor = DAPI (fluorescent DNA dye).

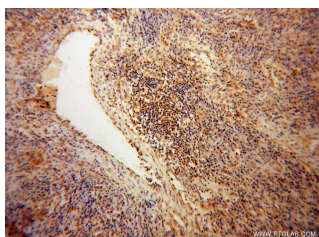

Immunohistochemical of paraffin-embedded human cervical cancer using 14307-1-AP (CDC16 antibody) at dilution of 1:50 (under 10x lens)

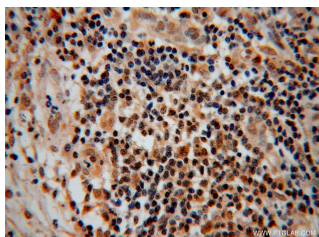

Immunohistochemical of paraffin-embedded human cervical cancer using 14307-1-AP (CDC16 antibody) at dilution of 1:50 (under 40x lens)

Applications

Tested applications:	ELISA, WB, IHC, IF
Species specificity:	Human, Mouse, Rat; other species not tested.
Calculated CDC16/APC6 MW:	72kd
Observed CDC16/APC6 MW:	72kd
Positive WB detected in	Jurkat cells, HeLa cells, K-562 cells
Positive IHC detected in	Human cervical cancer tissue
Positive IF detected in	HeLa cells
Recommended dilution:	WB: 1:200-1:2000 IHC: 1:20-1:200 IF: 1:20-1:200

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

Immunogen:	Ag5427
GenBank accession number:	BC010875
Gene ID (NCBI):	8881
Full name:	Cell division cycle 16 homolog (S. cerevisiae)

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3. Store at -20°C.