

ABCF1 Polyclonal Antibody

Catalog number: 13950-1-AP

Size: 30 µg/150 µl


Source: Rabbit

Isotype: IgG


Synonyms:

ABCF1; ABC27, ABC50, ABCF1,


ATP binding cassette 50

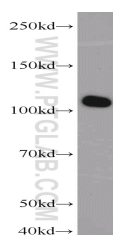

K-562 cells were subjected to SDS PAGE followed by western blot with 13950-1-AP(ABCF1 antibody) at dilution of 1:1500

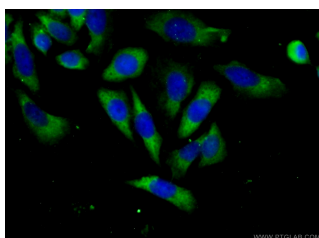

Immunofluorescent analysis of (-20% Ethanol) fixed HeLa cells using 13950-1-AP(ABCF1 Antibody) at dilution of 1:50 and Alexa Fluor 488-conjugated AffiniPure Goat Anti-Rabbit IgG(H+L)

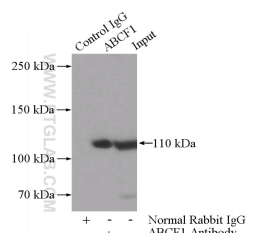

IP Result of anti-ABCF1 (IP:13950-1-AP, 4µg; Detection:13950-1-AP 1:500) with K-562 cells lysate 1200µg.

Applications

Tested applications:	ELISA, WB, IP, IF
Species specificity:	Human,Mouse,Rat; other species not tested.
Calculated ABCF1 MW:	92 kDa
Observed ABCF1 MW:	110 kDa
Positive WB detected in	K-562 cells
Positive IP detected in	K-562 cells
Positive IF detected in	HeLa cells
Recommended dilution:	WB: 1:500-1:5000
	IP: 1:200-1:2000
	IF: 1:20-1:200

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

Immunogen:	Ag5061
GenBank accession number:	BC034488
Gene ID (NCBI):	23
Full name:	ATP-binding cassette, sub-family F (GCN20), member 1

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3. Store at -20°C.