

## ABT1 Polyclonal Antibody

Catalog number: 14148-1-AP


Size: 71 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:

ABT1; ABT1, hABT1

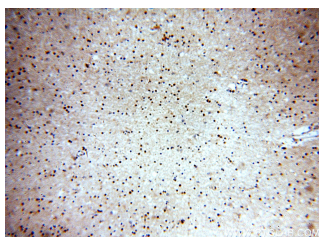

Immunohistochemical of paraffin-embedded human brain using 14148-1-AP (ABT1 antibody) at dilution of 1:100 (under 10x lens)

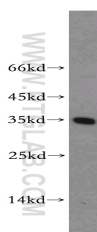

mouse liver tissue were subjected to SDS PAGE followed by western blot with 14148-1-AP (ABT1 antibody) at dilution of 1:300

### Applications

|                          |                                              |
|--------------------------|----------------------------------------------|
| Tested applications:     | ELISA, WB, IHC                               |
| Species specificity:     | Human, Mouse, Rat; other species not tested. |
| Calculated ABT1 MW:      | 31 kDa                                       |
| Observed ABT1 MW:        | 31-32 kDa                                    |
| Positive WB detected in  | Mouse liver tissue, HeLa cells               |
| Positive IHC detected in | Human brain tissue                           |
| Recommended dilution:    | WB: 1:200-1:1000<br>IHC: 1:20-1:200          |

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

### Immunogen information

|                           |                                    |
|---------------------------|------------------------------------|
| Immunogen:                | Ag5317                             |
| GenBank accession number: | BC066313                           |
| Gene ID (NCBI):           | 29777                              |
| Full name:                | Activator of basal transcription 1 |

### Product information

|                      |                                                                      |
|----------------------|----------------------------------------------------------------------|
| Purification method: | Antigen affinity purification                                        |
| Storage:             | PBS with 0.02% sodium azide and 50% glycerol pH 7.3. Store at -20°C. |