

CHST9 Polyclonal Antibody

Catalog number: 13135-1-AP

Size: 24 µg/150 µl


Source: Rabbit

Isotype: IgG

Synonyms:

CHST9; FLJ42328; GALNAC4ST-

2

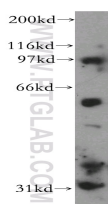

HepG2 cells were subjected to SDS PAGE followed by western blot with 13135-1-AP(CHST9 Antibody) at dilution of 1:1000

Applications

Tested applications:	ELISA, WB
Species specificity:	Human, mouse, rat; other species not tested.
Positive WB detected in:	HepG2 cells
Calculated CHST9 MW:	438aa, 52kd
Observed CHST9 MW:	52kd and 42kd
Recommended dilution:	WB: 1:500-1:5000

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP =

Immunoprecipitation, FC = Flow cytometry

Immunogen information

Immunogen:	Ag3874
GenBank accession number:	BC025764
Gene ID (NCBI):	83539
Full name:	Carbohydrate (N-acetylgalactosamine 4-O) sulfotransferase 9

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3. Store at -20°C.