

## Prepare Solutions with Confidence

### Simple and Affordable



#### Easy Dispensing Applications

The instrument automatically prepares your standards. Simply input the weight of your sample and set a target concentration. Minimize manual transcription errors by connecting a MS, ML, or ME balance for automatic weight transfer.



#### Clear Guidance

Designed for simplicity, the Easy Dose starts your routine task with a tap on the home screen. The user interface speaks your language – and 14 others. Learning and using a new analytical instrument has never been easier.




#### Simple Operation

Operate the instrument with the touchscreen, or the new responsive USB wheel controller. The wheel synchronizes effortlessly with the 4.3 inch touchscreen, offering exact control and a comfortable operation.



#### iTitrate for Easy Programming

iTitrate Intelligence makes method programming simple with EasySetup guidance. Set only the minimal parameters and the instrument programs the rest, including results calculated with your desired units.

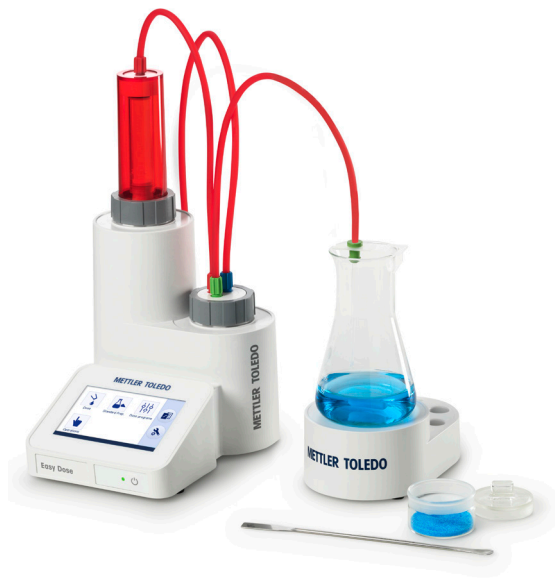

#### Easy Dose

Simplify and speeds up your everyday standard preparation and dispensing tasks. Dispensing tasks can be based on volume, speed, and time to tailor applications depending on your requirements. The simple user interface increases efficiency and simplicity. iTitrate Intelligence's guided set-up assists in method programming to ensure samples are prepared consistently, regardless of operator, and with minimal effort. Automatically export your results to a USB flash drive, or print them directly using a USB printer.

## Easy Dose

### Standard Delivery

Product	Order Number
EasyPlus burette 50 ml	30043902
EasyPlus tubing set instrument	30065464
EasyPlus drying tube	30044701
Valve	30042860
Bottle head	30060023
Burette cover set	30057633

### Accessories

Product	Order Number
EasyPlus burette 10 mL	30043900
EasyPlus burette 20 mL	30043901
EasyStir GT	30065467
EasyPlus titration head GT	30041102
Insert set titration head GT	30057635
Bottle adaptor China	30060026
Bottle adaptor India	30060024
Bottle adaptor Merck	23774
Bottle adaptor Fisher	23787
Siphon tips (5 pcs.)	23240
Glass beakers (80 mL), 20 pcs.	101446
Plastic beakers (100 mL), 120 pcs.	51109388
USB-P25 stripe printer	11124301
EasyWheel	30459699
UMix	30040000
UPlace	30019823
Tube holder, PP (3 pcs.)	30459710

### Technical Specifications

Power supply	Input voltage	24 VDC / 1.25 A
	Line input voltage	100 – 240 V $\pm$ 10%
	Input frequency	50 / 60 Hz
Environment	Ambient temperature	5 °C – 40 °C
	Atmospheric humidity	Max. 80% at 31 °C
Dimensions	Width x depth x height / weight	170 x 220 x 350 mm / 1.9 kg
Control and display	Touchscreen	4.3 inch, 480 x 272 pixel RGB
Communication	USB-A	USB full / low speed
	USB-B	USB full / low speed

[www.mt.com/EasyPlusTitration](http://www.mt.com/EasyPlusTitration)

For more information



#### METTLER TOLEDO Group

Local contact: [www.mt.com/contacts](http://www.mt.com/contacts)

Subject to technical changes  
© 06/2018 METTLER-TOLEDO. All rights reserved  
30469390  
Group MarCom 2479 LK