

ACBD3 Polyclonal Antibody

Catalog number: 14096-1-AP

Size: 26 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:


ACBD3; ACBD3, GCP60,
GOCAP1, Golgi phosphoprotein

1, Golgi resident protein


GCP60, GOLPH1, PAP7

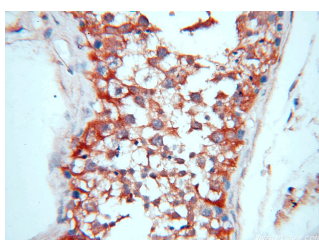


Immunohistochemical of paraffin-embedded human testis using 14096-1-AP (ACBD3 antibody) at dilution of 1:100 (under 40x lens)

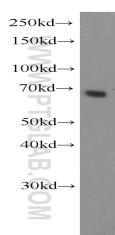

HeLa cells were subjected to SDS PAGE followed by western blot with 14096-1-AP (ACBD3 antibody) at dilution of 1:1200

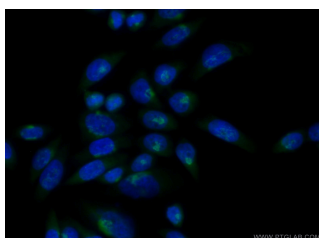

Immunofluorescent analysis of (10% Formaldehyde) fixed HeLa cells using 14096-1-AP (ACBD3 Antibody) at dilution of 1:50 and Alexa Fluor 488-conjugated AffiniPure Goat Anti-Rabbit IgG(H+L)

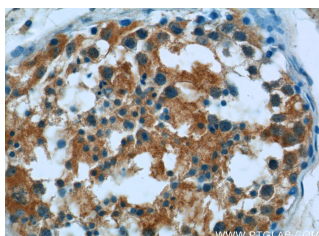

Applications

Tested applications:

ELISA, WB, IHC, IF

Species specificity:

Human, Mouse, Rat; other species not tested.

Calculated ACBD3 MW:

528aa, 61 kDa

Observed ACBD3 MW:

65-70 kDa

Positive WB detected in

HeLa cells, A2780 cells, Jurkat cells, K-562 cells, mouse ovary tissue, mouse testis tissue

Positive IHC detected in

Human testis tissue, human brain tissue, human heart tissue, human kidney tissue, human ovary tissue, human placenta tissue, human skin tissue, human spleen tissue

Positive IF detected in

HeLa cells

Recommended dilution:

WB: 1:500-1:5000

IHC: 1:20-1:200

IF: 1:50-1:500

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

Immunogen:

Ag5228

GenBank accession number:

BC060792

Gene ID (NCBI):

64746

Full name:

Acyl-Coenzyme A binding domain containing 3

Product information

Purification method:

Antigen affinity purification

Storage:

PBS with 0.02% sodium azide and 50% glycerol pH 7.3. Store at -20°C.

Immunohistochemistry of
paraffin-embedded human
testis tissue slide using 14096-
1-AP(ACBD3 Antibody) at
dilution of 1:50 (under 40x
lens)