

S. aureus LukS-PV (tag-free)

Catalog #: 0540-001

Lot #: 1503008

Description: Purified, tag-free *Staphylococcus aureus* Pantone-Valentine Leukocidin (PVL) S subunit (LukS-PV). The theoretical molecular weight of the protein is 32,465 Daltons.

Storage: 2-3 weeks at -20°C, -80°C long term

Size: 100 µg of protein is supplied in PBS at a concentration of 1.709 mg/mL. Protein demonstrates a molecular weight of approximately 35 kDa.

Relevance: This protein may be used in functional PVL toxicity assays in combination with LukF-PV, or as a control protein in ELISA assays or Western blotting when detecting LukS in PVL (+) strains of *S. aureus*.

Recommended Dilutions:

ELISA: Assay-dependent dilution.


WB: Assay-dependent dilution; internal QC demonstrates detection of 100 ng of LukS-PV protein when detected with anti-PVL LukS polyclonal antibody (cat# 04-0009) in Western blotting.

PVL Cytotoxicity assay: Cytotoxicity can be detected in human neutrophils when used in combination with LukF-PV in a concentration range of 1-100 nM.


For additional *S. aureus* products, please visit:

<http://ibtbioservices.com/index.php/product-and-reagents/staphylococcal-products>

SDS-PAGE and Western Blot Detection

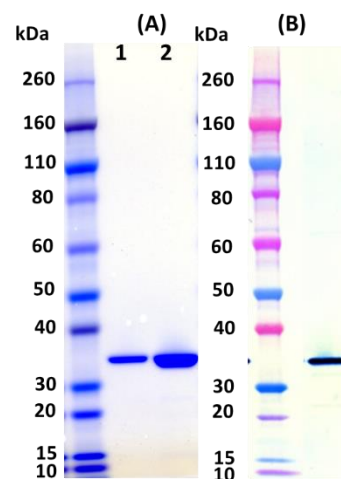

A) SDS-PAGE of LukS-PV: 1 µg (lane 1) and 5 µg (lane 2). (B) Western blot detection of LukS-PV at 100 ng, using IBT's anti-LukS-PV polyclonal antibody (Cat#04-0009) at 0.5 µg/mL and an anti-rabbit IgG-HRP conjugate, followed by TMB substrate.

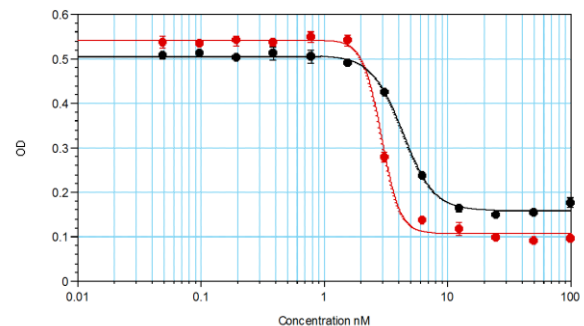

Toxin Functionality: Human promyelocytic leukemia cell line HL60 was differentiated into neutrophils by treatment with DMSO. Neutrophils were incubated with serial dilutions of LukS-PV and LukF-PV at equimolar concentration for 3 hours at 37°C with 5% CO₂ and 95% humidity. Cellular viability was determined by adding XTT and incubating for additional 16 hours. Cells were centrifuged and the OD determined in the supernatants at 470/690 nm. EC₅₀ values were found to be 2.92 nM for the LukS-PV tag-free lot 1503008 (red circles) and 4.38 nM for the LukS-PV his-tag lot 1108003 (black circles).

Intended for research use only. Not for human, therapeutic, or diagnostic applications.

The buyer cannot sell or otherwise transfer this product for Commercial Purposes without written approval of Integrated BioTherapeutics, Inc.

Copyright 2015. Integrated BioTherapeutics, Inc. All rights reserved