

Diaphragm Pumps for Air, Gases and Vapors

INNOVATIVE
TECHNOLOGY
WORLDWIDE



NEUBERGER, INC.



LABOPORT® N842.3FTP Pump

LABOPORT® Diaphragm Vacuum Pump

Technical features:

- 100% oil-free transfer
- Pure transfer, evacuation and compression
- Highly compatible with vapors and condensation
- Chemically-resistant
- Therefore suitable for highly aggressive or corrosive gases and vapors
- Maintenance-free
- Environmentally friendly
- Gastight, leakage rate approx. 6×10^{-3} mbar x l/s, not tested in serial production.

The chemically-resistant series N842.3 diaphragm pump is a double-head, oil-free device used in a wide range of laboratory applications. They transfer and pump down without contamination.

The heart of this very compact pump is a KNF structured diaphragm. This patented diaphragm was stress-optimized using the Finite Elements method. As a result, we were able to make the pump smaller while increasing the service life of the diaphragm.

Technical data:

N842.3FTP

Delivery (l/min) ¹⁾	34
Ultimate vacuum (Torr)	1.5
Operating pressure (psi)	15
Connectors for tube (in.)	ID 3/8"
Permissible gas and ambient temperature	+5...+40 °C
Voltage/Frequencies	115V/60Hz
Motor protection	IP 44
Power P ₁	290 W
Operating current	3.2 A
Weight	27 lbs.
Dimensions	
LxHxW (mm)	341/223/167
With thermal switch and power fuse	

Motors with other voltages and frequencies on request.

Material in contact with the pumped media

Type/OrderNo.	Pump head	Diaphragm	Valves
N842.3FTP	PTFE	PTFE-coated	FFPM

¹⁾ at atm. pressure

Diaphragm Pumps for Air, Gases and Vapors


INNOVATIVE
TECHNOLOGY
WORLDWIDE




NEUBERGER, INC.

Dimensions and performance characteristics

Pump down time for 10 l receiver

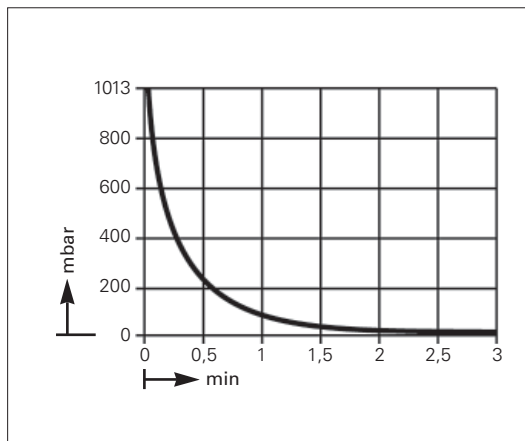

Performance characteristics

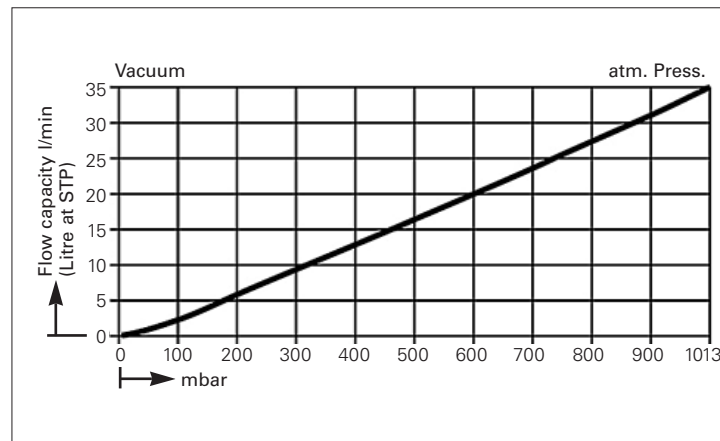

Dimensions (mm)

