


Description: Microporous sleeve covers (Model No. XC6618)
Size: 18"
Packaging: 20/PK; 10 PKS/CS; 200/CS
Color: White
Material: Microporous 60g, antistatic
 Latex-free elastic

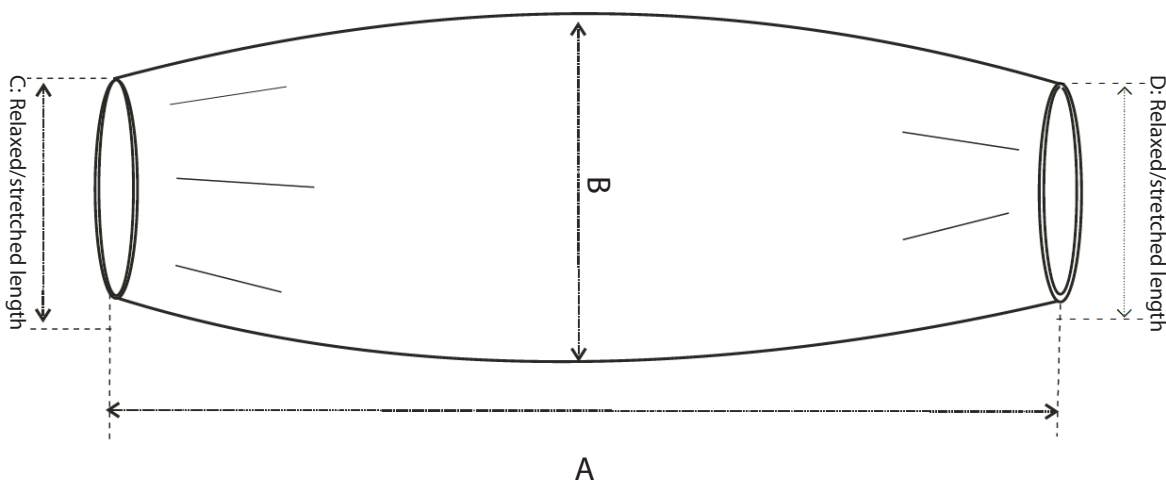

Item	18"	Tolerance	Comment
A	46cm	+/-2cm	
B	20cm	+/-2cm	Width
C	8/20cm	+/-2cm	Relaxed/stretched length
D	7.5/19.5cm	+/-2cm	Relaxed/stretched length

Material Inspection Report

XC6000 Material

Physical Properties

Test Item	Test Standard	Result
Weight (g/m²)	ASTM D3776/D-3776M-20 (Option C))	57.7
Ball Burst Strength, lbf	ASTM D3787-16	105.0
Mullen Burst Strength, psi	ASTM D3786/D3786M-18	32
Tearing Strength Length, lbf Width, lbf	ASTM D5587-15	11.5 4.8
Tensile Strength Length, lbf Width, lbf	ASTM D 5034-09, Grab Method	29.2 18.4
Water Vapor Transmission, g/24h/m²	ASTM E96/E96M-16 (Procedure B)	878
Particle Shedding (Helmke Drum)	IEST-RP-CC003.4	Category II
Electrostatic Decay of Fabric	AATCC TM 76-2000e3(2018)e	PASS

Barrier Properties

Water Resistance, Impact	AATCC 42-2017e	0.00
Synthetic Blood Resistance Penetration, 100%	ASTM F1670/F1670M-2017, (Procedure B)	PASS

Chemical Resistance

Household Bleach 100%	ASTM F903-03	PASS
Hydrogen Peroxide	ASTM F903-03	PASS
Phosphoric Acid 85%	ASTM F903-03	PASS
Sulfuric Acid 70%	ASTM F903-03	FAIL