

BLK Polyclonal Antibody

Catalog number: 10510-1-AP


Size: 36 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:


BLK; B lymphocyte kinase, B lymphoid tyrosine kinase, BLK, p55 Blk, Tyrosine protein kinase Blk

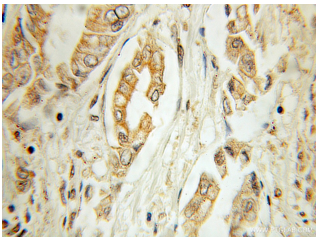

Immunohistochemical of paraffin-embedded human pancreas cancer using 10510-1-AP (BLK antibody) at dilution of 1:10 (under 10x lens)

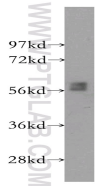

SH-SY5Y cells were subjected to SDS PAGE followed by western blot with 10510-1-AP (BLK antibody) at dilution of 1:500

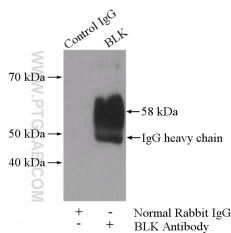

IP Result of anti-BLK (IP:10510-1-AP, 4µg; Detection:10510-1-AP 1:500) with SH-SY5Y cells lysate 1600µg.

Applications

Tested applications:	ELISA, WB, IHC, IP
Species specificity:	Human, Mouse, Rat; other species not tested.
Calculated BLK MW:	58 kDa
Observed BLK MW:	58 kDa
Positive WB detected in:	SH-SY5Y cells
Positive IP detected in:	SH-SY5Y cells
Positive IHC detected in:	Human pancreas cancer tissue
Recommended dilution:	WB: 1:500-1:5000
	IP: 1:200-1:2000
	IHC: 1:10-1:100

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

Immunogen:	Ag0795
GenBank accession number:	BC007371
Gene ID (NCBI):	640
Full name:	B lymphoid tyrosine kinase

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C.