

DENR Polyclonal Antibody

Catalog number: 10656-1-AP

Size: 22 µg/150 µl


Source: Rabbit

Isotype: IgG


Synonyms:

DENR; DENR, density regulated protein, DRP, DRP1, Protein

DRP1, SMAP 3

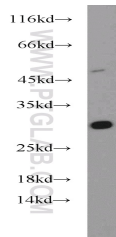

A431 cells were subjected to SDS PAGE followed by western blot with 10656-1-AP (DENR antibody) at dilution of 1:1000

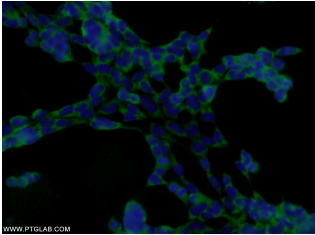

Immunofluorescent analysis of HEK-293 cells using 10656-1-AP (DENR Antibody) at dilution of 1:25 and Alexa Fluor 488-conjugated AffiniPure Goat Anti-Rabbit IgG(H+L)

Applications

Tested applications:	ELISA, WB, IF
Species specificity:	Human, Mouse, Rat; other species not tested.
Calculated DENR MW:	27 kDa
Observed DENR MW:	27-30kd
Positive WB detected in	A431 cells, HEK-293 cells, mouse brain tissue, mouse lung tissue
Positive IF detected in	HEK-293 cells, A431 cells
Recommended dilution:	WB: 1:500-1:5000 IF: 1:10-1:100

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

Immunogen:	Ag1023
GenBank accession number:	BC007860
Gene ID (NCBI):	8562
Full name:	Density-regulated protein

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C.