

## NOL8 Polyclonal Antibody

Catalog number: 12043-1-AP

Size: 35 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:

NOL8; C9orf34;

DKFZp686P12242; FLJ20736;

Nop132; bA62C3.3; bA62C3.4

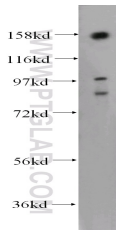

SGC-7901 cells were subjected to SDS PAGE followed by western blot with 12043-1-AP (NOL8 Antibody) at dilution of 1:1200



Immunofluorescent analysis of HepG2 cells, using NOL8 antibody 12043-1-AP at 1:25 dilution and Rhodamine-labeled goat anti-rabbit IgG (red). Blue pseudocolor = DAPI (fluorescent DNA dye).

### Applications

Tested applications:	ELISA, WB, IF
Species specificity:	Human, mouse, rat; other species not tested.
Positive WB detected in:	SGC-7901 cells
Calculated NOL8 MW:	132kd
Observed NOL8 MW:	150kd
Positive IF detected in:	HepG2 cells
Recommended dilution:	WB: 1:500-1:5000 IF: 1:20-1:50

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation, FC = Flow cytometry

### Immunogen information

Immunogen:	Ag2672
GenBank accession number:	BC013788
Gene ID (NCBI):	55035
Full name:	Nucleolar protein 8

### Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3. Store at -20°C.