

RAB40B Polyclonal Antibody

Catalog number: 12471-1-AP


Size: 42 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:

RAB40B; FLJ42385, Protein Rar, RAB40B, RAR, Ras related protein Rab 40B, SEC4L

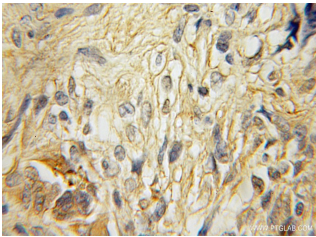

Immunohistochemical of paraffin-embedded human gliomas using 12471-1-AP (RAB40B antibody) at dilution of 1:50 (under 10x lens)

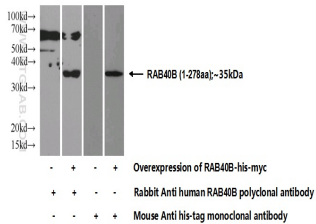

Transfected HEK-293 cells were subjected to SDS PAGE followed by western blot with 12471-1-AP (RAB40B Antibody) at dilution of 1:1000

Applications

Tested applications:	ELISA, IHC, WB
Species specificity:	Human; other species not tested.
Calculated RAB40B MW:	278aa, 31 kDa
Observed RAB40B MW:	
Positive WB detected in	Transfected HEK-293 cells
Positive IHC detected in	Human gliomas tissue
Recommended dilution:	WB: 1:500-1:5000 IHC: 1:20-1:200

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

Immunogen:	Ag3155
GenBank accession number:	BC018039
Gene ID (NCBI):	10966
Full name:	RAB40B, member RAS oncogene family

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C.