



Coriolis Biological Air Samplers

Air bio-contamination control


- Indoor & Outdoor sampling
- Collect any biological sample
- Compatible with multiple downstream analysis

Air contamination: a worldwide issue for public health

The Coriolis air samplers can collect biological particles in the air, offering new perspectives for the monitoring of airborne contamination.

- A versatile range of instruments to meet your needs
- Patented technologies adapted to the collection of viruses, bacteria, molds, pollens, spores
- Sample output compatible with standard methods for culture and molecular biology

WORKFLOW

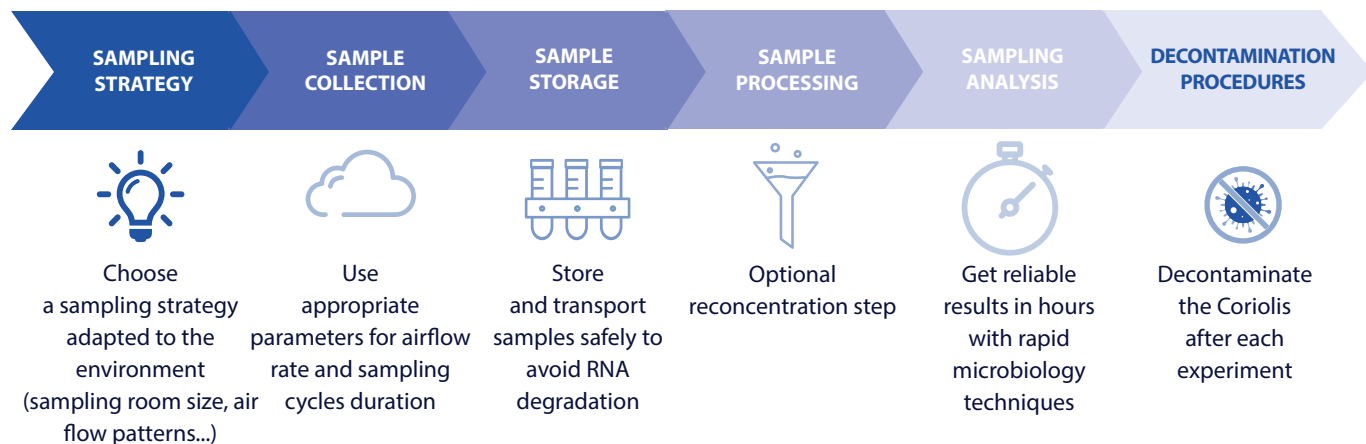


RESEARCH AREAS



Pollution & Environment



Biomedical & Health

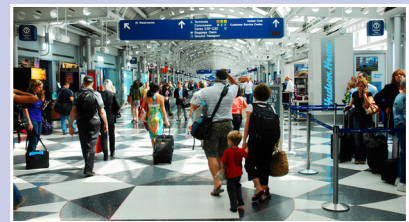

Indoor Air Quality



Food, Pharma, Veterinary & Industry



Research & Development



Defense & Civil Security

Coriolis:


a range of solutions for air monitoring

The Coriolis range was developed to go beyond the traditional technical performances of airborne bio-contamination control.



CORIOLIS+

AIR SAMPLER COMPATIBLE
WITH ANY TYPE OF ANALYSIS

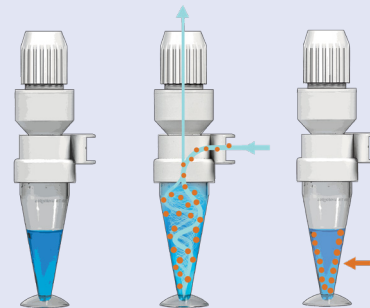


- The cone is prefilled with a collection liquid
- Aspirated air runs into the cone and forms a vortex
- Particles are centrifuged into the cone walls and trapped into the collection liquid
- Particles in the liquid sample are ready for analysis



CORIOLIS COMPACT

PORTABLE DRY AIR SAMPLER
FOR BIOCONTAMINATION CONTROL

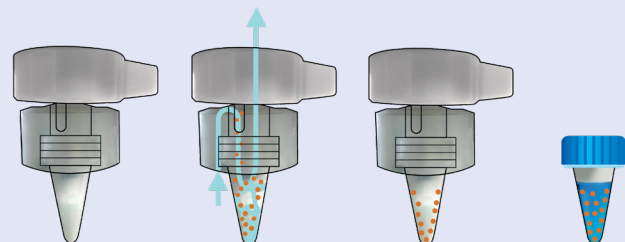

- The cone is placed into the device and locked
- Aspirated air runs into the cone and forms a vortex
- Particles in the air are centrifuged into the cone walls
- Collected particles are recovered by rinsing the cone
- Remote control with Bertin Connect Application

COROLIS+

P0010080-CORM0-B

Technology: wet

Dimensions: 220 x 330 x 360 mm
Weight: 2.8 kg.

Noise: 66dBa at 300L/min

Use temperature: between + 5 / + 40°C

Autonomy: 2 h
Battery: Li-Ion

Air flow: 100-300 L/min

Collect any biological sample

Use: indoor & outdoor

Downstream analysis: multiple
> Culture/PCR



COROLIS COMPACT

P002055-CORC0-A

Technology: dry

Dimensions: 255 x 135 x 130 mm
Weight: 1.42 kg.

Noise: 62dBa

Use temperature: between +5 / +45°C

Autonomy: max 8 h at 20°C
Battery: Li-Ion

Air flow: 50 L/min

Collect any biological sample


Use: indoor & outdoor
Bluetooth

Downstream analysis: multiple
> Culture/PCR



Discover our solutions

Bertin Life Sciences and Bertin Bioreagent: Comprehensive Life Sciences Solutions

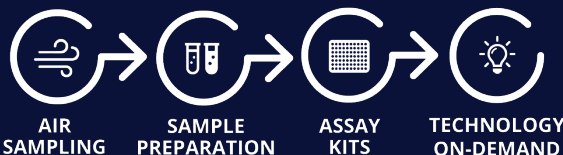

HEAD OFFICE
Bertin Technologies SAS
Parc d'activités du Pas du Lac
10 bis, avenue Ampère
78 180 Montigny-le Bretonneux
FRANCE

US OFFICE
Bertin Corp.
11630-11634 Boiling Brook Parkway,
Rockville MD 20852, United States

ORDER
+33 (0)1 39 30 60 36
order-life@bertin.group
order@bertin-corp.com (US)

TECHNICAL SUPPORT & MAINTENANCE
maintenance-life@bertin.group
air-sampler@bertin.group