

SNX4 Polyclonal Antibody

Catalog number: 11506-1-AP


Size: 35 µg/150 µl

Source: Rabbit

Isotype: IgG

Synonyms:

SNX4;

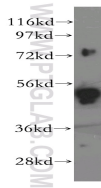

NIH/3T3 cells were subjected to SDS PAGE followed by western blot with 11506-1-AP(SNX4 Antibody) at dilution of 1:1000

Applications

Tested applications:	ELISA, WB
Species specificity:	Human, mouse, rat; other species not tested.
Positive WB detected in:	NIH/3T3 cells
Calculated SNX4 MW:	450aa,52kd
Observed SNX4 MW:	50kd
Recommended dilution:	WB: 1:500-1:5000

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation, FC = Flow cytometry

Immunogen information

Immunogen:	Ag2058
GenBank accession number:	BC018762
Gene ID (NCBI):	8723
Full name:	Sorting nexin 4

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C.