

ST3GAL2 Polyclonal Antibody

Catalog number: 13620-1-AP

Size: 21 µg/150 µl

Source: Rabbit

Isotype: IgG


Synonyms:

ST3GAL2; Alpha 2,3 ST 2, Gal

NAc6S, Sialyltransferase 4B,

SIAT4 B, SIAT4B, ST3Gal II,

ST3GAL2, ST3GalA.2, ST3GALII

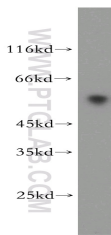

human heart tissue were subjected to SDS PAGE followed by western blot with 13620-1-AP(ST3GAL2 antibody) at dilution of 1:500

Applications

Tested applications:	ELISA, WB
Species specificity:	Human,Mouse,Rat; other species not tested.
Calculated ST3GAL2 MW:	350aa,40 kDa
Observed ST3GAL2 MW:	55kd
Positive WB detected in	Human heart tissue, HL-60 cells, human colon tissue, human liver tissue, human spleen tissue, K-562 cells
Recommended dilution:	WB: 1:500-1:5000

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

Immunogen:	Ag4528
GenBank accession number:	BC036777
Gene ID (NCBI):	6483
Full name:	ST3 beta-galactoside alpha-2,3-sialyltransferase 2

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3. Store at -20°C.