

AGR3 Polyclonal Antibody

Catalog number: 11967-1-AP


Size: 30 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:


AGR3; AG 3, AG3, AGR3, hAG 3, HAG3, PDIA18

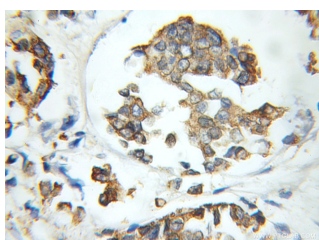

Immunohistochemical of paraffin-embedded human breast cancer using 11967-1-AP (AGR3 antibody) at dilution of 1:50 (under 10x lens)

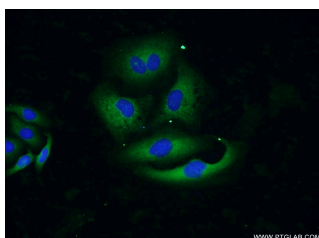

Immunofluorescent analysis of SKOV-3 cells using 11967-1-AP (AGR3 Antibody) at dilution of 1:50 and Alexa Fluor 488-conjugated AffiniPure Goat Anti-Rabbit IgG(H+L)

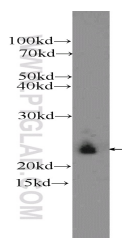

SKOV-3 cells were subjected to SDS PAGE followed by western blot with 11967-1-AP (AGR3 Antibody) at dilution of 1:600

Applications

Tested applications:	ELISA, IHC, WB, IF
Species specificity:	Human, Mouse; other species not tested.
Calculated AGR3 MW:	166aa, 19 kDa
Observed AGR3 MW:	19-25 kDa
Positive WB detected in	SKOV-3 cells, MCF-7 cells
Positive IHC detected in	Human breast cancer tissue
Positive IF detected in	SKOV-3 cells
Recommended dilution:	WB: 1:200-1:2000 IHC: 1:20-1:200 IF: 1:20-1:200

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

Immunogen:	Ag2577
GenBank accession number:	BC058284
Gene ID (NCBI):	155465
Full name:	Anterior gradient homolog 3 (Xenopus laevis)

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C.