

SHKBP1 Polyclonal Antibody

Catalog number: 12703-1-AP


Size: 23 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:

SHKBP1; PP203; Sb1

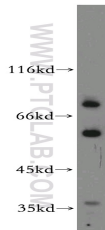

A431 cells were subjected to SDS PAGE followed by western blot with 12703-1-AP(SHKBP1 Antibody) at dilution of 1:600

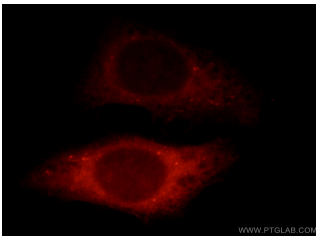

Immunofluorescent analysis of HepG2 cells, using SHKBP1 antibody 12703-1-AP at 1:25 dilution and Rhodamine-labeled goat anti-rabbit IgG (red).

Applications

Tested applications:	ELISA, WB, IF
Species specificity:	Human, mouse, rat; other species not tested.
Positive WB detected in:	A431 cells
Calculated SHKBP1 MW:	707aa,76kd
Observed SHKBP1 MW:	76kd
Positive IF detected in:	HepG2 cells
Recommended dilution:	WB: 1:200-1:2000 IF: 1:20-1:50

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation, FC = Flow cytometry

Immunogen information

Immunogen:	Ag3386
GenBank accession number:	BC022855
Gene ID (NCBI):	92799
Full name:	SH3KBP1 binding protein 1

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3. Store at -20°C.