

FMR1NB Polyclonal Antibody

Catalog number: 11069-2-AP


Size: 22 µg/150 µl

Source: Rabbit

Isotype: IgG

Synonyms:

FMR1NB; CT37; FLJ25736; NY-SAR-35

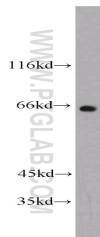

MCF7 cells were subjected to SDS PAGE followed by western blot with 11069-2-AP(FMR1NB Antibody) at dilution of 1:800

Applications

Tested applications:	ELISA, WB
Species specificity:	Human; other species not tested.
Positive WB detected in:	MCF7 cells, A549 cells, SKOV-3 cells
Calculated FMR1NB MW:	29kd
Observed FMR1NB MW:	66kd
Recommended dilution:	WB: 1:200-1:2000

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation, FC = Flow cytometry

Immunogen information

Immunogen:	Ag1531
GenBank accession number:	BC034320
Gene ID (NCBI):	158521
Full name:	Fragile X mental retardation 1 neighbor

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C.