

Sieve Cleaning Module SRH 4/200 with integrated turning device




Elma Order Numbers

Sieve cleaning module 230 V	104 2706
Sieve cleaning module 115 V	104 4904

Sieve cleaning module SRH 4/200 with integrated turning device for Elma ultrasonic units


- Time and cost saving intensive cleaning of sieves
- Cleaning of up to 4 sieves [Ø max. 200/203 mm (8")] in one cleaning cycle
- Sieves of different brands can be cleaned simultaneously
- Easy handling, the module is placed into the ultrasonic tank like a basket
- Metal components made of stainless steel, hygienic and easy to clean
- Inclined position of the sieves guarantees improved cleaning and rinsing results
- Removed dirt particles collect on the tank floor, not in the sieves
- Gentle ultrasonic cleaning; mesh size, grid tension or measuring accuracy are not changed

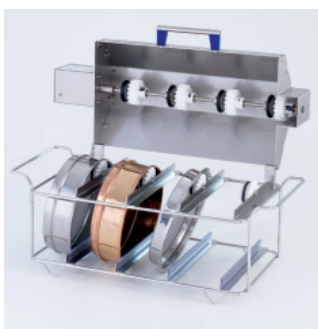

Sieve cleaning module loaded



Sieve cleaning module inserted in ultrasonic unit

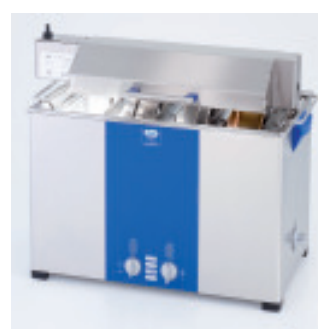

Sieve cleaning module operating



Sieve cleaning module placed onto plastic cover (Elmasonic 300) for drip-off

Technical Details

Mains voltage (Vac)	230 or 115
Mains frequency (Hz)	50 / 60
Outer dimensions W / D / H (mm)	580 / 260 / 300
Weight (ca. kg)	6.8
Material	stainless steel

Suitable for Elma ultrasonic units

Elmasonic S 300 (H) / P 300 H	for 1 SRH 4/200
Elmasonic S 450 H	for 1 SRH 4/200
Elmasonic S 900 H	for 2 SRH 4/200