

labForce[®] Amber Sterile HDPE Wide Mouth Bottles

labForce[®] Sterile Amber HDPE Wide Mouth bottles are made of high-quality HDPE raw material, and produced in a Class 100,000 grade cleanroom production environment with multiple quality system certifications. Made from natural HDPE, which is non-biotoxic, these bottles have excellent physical and chemical properties, are highly resistant to leaks, liquid impact and acids, resistant to low temperatures of -40°C and meets the requirements of autoclave processing. Amber bottles are resistant to organic solvents and ideal for storing photosensitive substances.



HIGH QUALITY CONSTRUCTION

Blow-molded with no seams, the advanced molding technology ensures that the bottles have a more uniform wall thickness, resulting in high quality walls to prevent breakage or puncturing of the bottom, and a flat bottom for smooth running on the filling line.

CLEAN AND STERILE

STERILE labForce[®] sterile HDPE bottles are free of DNase, RNase, protease, and exogenous DNA/RNA, and are non-pyrogenic. Sterilized by electron beam irradiation.

LEAK-PROOF DESIGN

labForce[®] bottle closures are designed with a liquid-proof mouth so that repeated opening and tightening will not affect the sealing performance of the bottle. The sealing ring inside the closure fits snugly against the beveled inner edge of the bottle neck to form a leak-proof valve to ensure the closure is sealed. labForce[®] bottle closures, made of polypropylene are designed with a liquid-proof mouth so that repeated opening and tightening will not affect the sealing performance of the bottle. The sealing ring inside the closure fits snugly against the beveled inner edge of the bottle neck to form a leak-proof valve to ensure the closure is sealed.

COMPLETE RANGE OF SIZES


30ml, 60ml, 125ml, 250ml, 500ml, 1000,ml

Thomas No.	Description	Size	Packs	Case
21A00K990	labForce [®] Amber Sterile HDPE Wide Mouth Bottles Closures Attached	30ml	10/packs; 20 packs/case	200
21A00K991	labForce [®] Amber Sterile HDPE Wide Mouth Bottles Closures Attached	60ml	10/packs; 10 packs/case	200
21A00K992	labForce [®] Amber Sterile HDPE Wide Mouth Bottles Closures Attached	125ml	10/packs; 10 packs/case	100
21A00K993	labForce [®] Amber Sterile HDPE Wide Mouth Bottles Closures Attached	250ml	10/packs; 10 packs/case	100
21A00K994	labForce [®] Amber Sterile HDPE Wide Mouth Bottles Closures Attached	500ml	5/packs; 10 packs/case	50
21A00K995	labForce [®] Amber Sterile HDPE Wide Mouth Bottles Closures Attached	1000ml	3/packs; 9 packs/case	27



Manufactured in ISO 13485:2016 Facility

Wide mouth bottle specification						
Content Name	Bottle Material	D1(mm)	H1(mm)	D2(mm)	D3(mm)	H2(mm)
30ml	PE	Ø 33.9	59.05	Ø 20.9	Ø 32.2	15.5
60ml	PE	Ø 38.6	81.65	Ø 20.9	Ø 32.2	15.5
125ml	PE	Ø 50.8	94.7	Ø 27.8	Ø 42.2	15.5
250ml	PE	Ø 61.7	127.6	Ø 33.0	Ø 46.4	16.8
500ml	PE	Ø 73.0	161.6	Ø 43.8	Ø 57.5	16.8
1000ml	PE	Ø 92.2	193.9	Ø 43.7	Ø 57.5	16.8

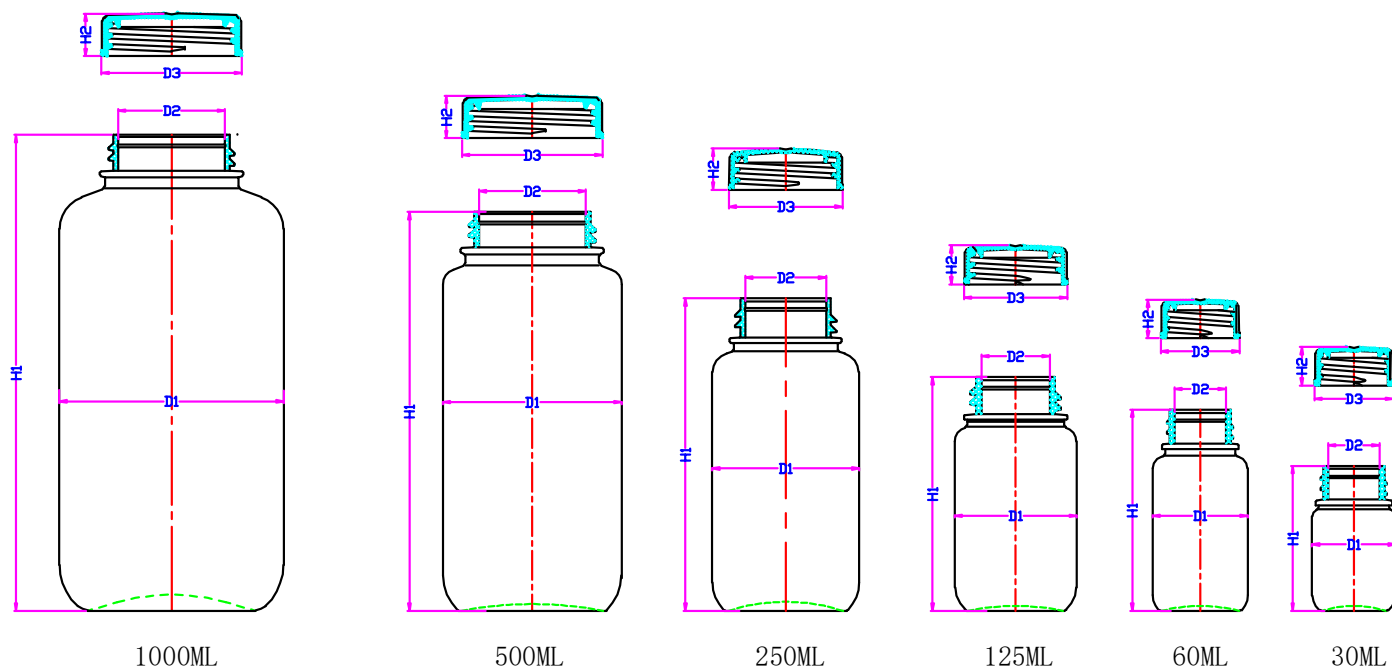

REV 4/11/22