




| Width                   | 23 7/8"   |  |
|-------------------------|---|--|
| Interior Finish         | White   |  |
| Total Capacity          | 5.8 Cubic Feet; 166 Liters  |  |
| Refrigerator Capacity   | 4.4 Cubic Feet; 124 Liters  |  |
| Freezer Capacity        | 1.4 Cubic Feet; 40 Liters   |  |
| Toe Grill Finish        | Black   |  |
| Hinge                   | Cam-assisted self close   |  |
| Handle                  | Discrete  |  |
| Lock                    | Yes   |  |
| Temperature Range       | Refrigerator: 34°F to 45°F (1°C to 7°C)<br>Freezer: 0°F to 15°F (-18°C to -9°C) |  |
| Electrical Requirements | 115V / 60Hz / 3.3A  |  |
| Length of Power Cord    | 5'  |  |
| Product Dimensions      | 23 <sup>7</sup> / <sub>8</sub> " W x 34" H x 24" D                              |  |
| Shipping Weight         | 120 lbs; 142 lbs with skid  |  |
| Finishes                | White   |  |
| Agency                  | UL Listed   |  |

| Ordering Details                                    |            |
|---|------------|
| Solid White Cabinet and Door with Lock, Right Hinge | MS24RFS2RW |

## **Rough In Dimensions**

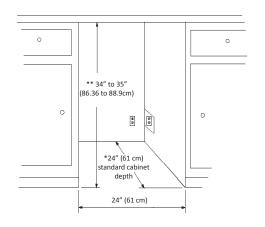

## 24" Refrigerator Freezer

Model # MS24RFS2RW

- Front-vented, permitting it to be built into casework or installed under lab benches without the need for additional ventilation clearance
- Knob-type electronic thermostat for selection of set point temperature is located in the toe-space grille below the door
- Accurate temperature control of contents within +3°/-3°C of set point.
- Recommended ambient temperature 18.3°C to 32.2°C (65°F to 90°F)
- Non-CFC 245FA blown foam insulation in cabinet walls and door, seals in the cold
- 22 gauge textured steel cabinet
- Solid door configurations have a "box style" textured steel exterior and three interior door shelves and a discrete top edge mounted door handle
- 4.4 cu.ft. of refrigerator space and 1.4 cu.ft. of freezer space
- Forced air condenser system efficiently removes heat and allows this unit to be built into casework without the loss of cooling performance
- Powerful compressor for fast pull-down time and temperature recovery after door opening or commodity loading
- Two adjustable and removable wire shelves
- "Box Style" door for easy mounting of inventory management system hardware
- Cam-assisted door close hinge system to reduce "door-ajar" events
- 4 Leveling legs for easy height adjustment
- Extra-long, jacketed commercial power cord with hospital grade plug to meet stringent hospital and clinical use standards
- Durable white powder coat paint finish standard
- Commercial UL Listed and warrantied for commercial use
- Probe port in cabinet wall for convenient central temperature monitoring probe installation
- Built-in door lock for secure storage

## **Product Dimensions**

