

PHACTR4 Polyclonal Antibody

Catalog number: 13408-1-AP

Size: 51 µg/150 µl


Source: Rabbit

Isotype: IgG


Synonyms:

PHACTR4; DKFZp686L07205,

FLJ13171, PHACTR4

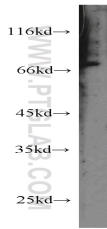

mouse eye tissue were subjected to SDS PAGE followed by western blot with 13408-1-AP(PHACTR4 antibody) at dilution of 1:400

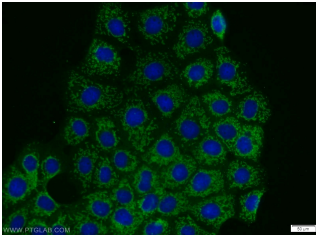

Immunofluorescent analysis of A431 cells using 13408-1-AP(PHACTR4 Antibody) at dilution of 1:50 and Alexa Fluor 488-conjugated AffiniPure Goat Anti-Rabbit IgG(H+L)

Applications

Tested applications:	ELISA, WB, IF
Species specificity:	Human,Mouse,Rat; other species not tested.
Calculated PHACTR4 MW:	702aa,78 kDa
Observed PHACTR4 MW:	71-85 kDa
Positive WB detected in	Mouse eye tissue, A431 cells, human liver tissue
Positive IF detected in	A431 cells
Recommended dilution:	WB: 1:200-1:2000 IF: 1:20-1:200

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

Immunogen:	Ag4054
GenBank accession number:	BC029266
Gene ID (NCBI):	65979
Full name:	Phosphatase and actin regulator 4

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3. Store at -20°C.