

UBP1 Polyclonal Antibody

Catalog number: 13856-1-AP

Size: 27 µg/150 µl


Source: Rabbit

Isotype: IgG

Synonyms:

UBP1; DKFZp686L1745; LBP-1B;

LBP-1a; LBP1A; LBP1B

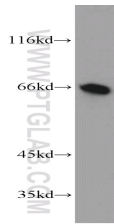

human kidney tissue were subjected to SDS PAGE followed by western blot with 13856-1-AP(UBP1 Antibody) at dilution of 1:1000

Applications

Tested applications:	ELISA, WB
Species specificity:	Human, mouse, rat; other species not tested.
Positive WB detected in:	Human kidney tissue
Calculated UBP1 MW:	60kd
Observed UBP1 MW:	68kd
Recommended dilution:	WB: 1:500-1:5000

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation, FC = Flow cytometry

Immunogen information

Immunogen:	Ag4830
GenBank accession number:	BC047235
Gene ID (NCBI):	7342
Full name:	Upstream binding protein 1 (LBP-1a)

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3. Store at -20°C.