

## PRAM1 Polyclonal Antibody

Catalog number: 13651-1-AP


Size: 22 µg/150 µl

Source: Rabbit

Isotype: IgG

Synonyms:

PRAM1; PML RAR, PRAM 1,  
PRAM1

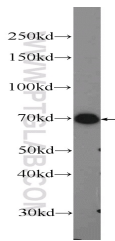

U-937 cells were subjected to SDS PAGE followed by western blot with 13651-1-AP (PRAM1 Antibody) at dilution of 1:300

### Background

PRAM1, also named as PML-RARA-regulated adapter molecule 1, is a 718 amino acid protein, which contains 1 SH3 domain. PRAM1 is expressed in peripheral blood leukocytes and bone marrow. It is not expressed in non hematopoietic tissues except in lung. PRAM1 may be involved in myeloid differentiation and integrin signaling in neutrophils.

### Applications

Tested applications:	ELISA, WB
Species specificity:	Human; other species not tested.
Calculated PRAM1 MW:	718aa, 79 kDa
Observed PRAM1 MW:	70 kDa
Positive WB detected in	U-937 cells
Recommended dilution:	WB: 1:200-1:1000

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

### Immunogen information

Immunogen:	Ag4162
GenBank accession number:	BC028012
Gene ID (NCBI):	84106
Full name:	PML-RARA regulated adaptor molecule 1

### Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C.