

APLF Polyclonal Antibody

Catalog number: 14252-1-AP

Size: 36 µg/150 µl


Source: Rabbit

Isotype: IgG

Synonyms:

APLF; APFL; C2orf13; FLJ16593;

MGC47799; PALF; Xip1

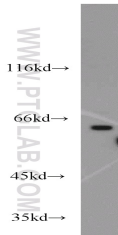

PC-3 cells were subjected to SDS PAGE followed by western blot with 14252-1-AP (APLF Antibody) at dilution of 1:600

Applications

Tested applications:	ELISA, WB
Species specificity:	Human; other species not tested.
Positive WB detected in:	PC-3 cells
Calculated APLF MW:	57kd
Observed APLF MW:	60kd
Recommended dilution:	WB: 1:200-1:2000

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation, FC = Flow cytometry

Immunogen information

Immunogen:	Ag5525
GenBank accession number:	BC041144
Gene ID (NCBI):	200558
Full name:	Aprataxin and PNKP like factor

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3. Store at -20°C.