

## Drebrin Polyclonal Antibody

Catalog number: 10260-1-AP


Size: 40 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:

DBN1; D0S117E, DBN1,  
DKFZp434D064, Drebrin,  
drebrin 1

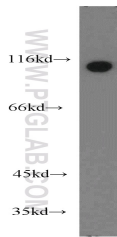

HeLa cells were subjected to SDS PAGE followed by western blot with 10260-1-AP (DBN1 antibody) at dilution of 1:1500

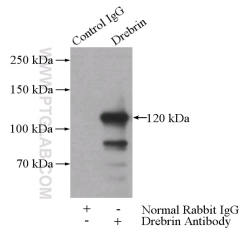

IP Result of anti-Drebrin (IP:10260-1-AP, 4µg; Detection:10260-1-AP 1:2000) with K-562 cells lysate 1200ug.

### Applications

Tested applications:	ELISA, WB, IP
Cited applications:	IF, IHC, WB
Species specificity:	Human, Mouse, Rat; other species not tested.
Cited species:	Mouse, rat
Calculated Drebrin MW:	71 kDa
Observed Drebrin MW:	100-130 kDa
Positive WB detected in:	HeLa cells, K-562 cells
Positive IP detected in:	K-562 cells
Recommended dilution:	WB: 1:500-1:5000 IP: 1:1000-1:10000

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

### Immunogen information

Immunogen:	Ag0371
GenBank accession number:	BC000283
Gene ID (NCBI):	1627
Full name:	Drebrin 1

### Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C.