

TES Polyclonal Antibody

Catalog number: 10258-1-AP


Size: 50 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:

TES; TES, Testin

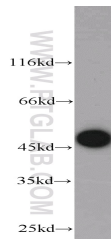

HEK-293 cells were subjected to SDS PAGE followed by western blot with 10258-1-AP(TES antibody) at dilution of 1:300

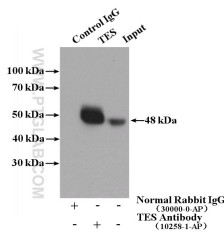

IP Result of anti-TES (IP:10258-1-AP, 4µg; Detection:10258-1-AP 1:1500) with K-562 cells lysate 4000µg.

Applications

Tested applications:	ELISA, WB, IP
Cited applications:	WB
Species specificity:	Human, Mouse, Rat; other species not tested.
Cited species:	Mouse
Calculated TES MW:	48 kDa
Observed TES MW:	38-48 kDa
Positive WB detected in	HEK-293 cells, K-562 cells, MCF7 cells, mouse testis tissue
Positive IP detected in	K-562 cells
Recommended dilution:	WB: 1:200-1:2000 IP: 1:500-1:5000

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

Immunogen:	Ag0379
GenBank accession number:	BC001451
Gene ID (NCBI):	26136
Full name:	Testis derived transcript (3 LIM domains)

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C.