

## CORO2B Polyclonal Antibody

Catalog number: 13802-1-AP


Size: 25 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:

CORO2B; CLIPINC; KIAA0925

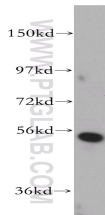

Jurkat cells were subjected to SDS PAGE followed by western blot with 13802-1-AP(CORO2B Antibody) at dilution of 1:1000

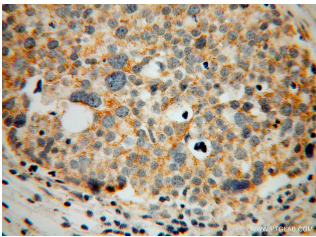

Immunohistochemistry of paraffin-embedded human ovary tumor using 13802-1-AP(CORO2B Antibody) at Dilution 1:100 (under 40x lens)

### Applications

Tested applications:	ELISA, WB, IHC
Species specificity:	Human, mouse, rat; other species not tested.
Positive WB detected in:	Jurkat cells, mouse brain tissue
Calculated CORO2B MW:	475aa,54kd
Observed CORO2B MW:	54kd
Positive IHC detected in:	Human ovary tumor
Recommended dilution:	WB: 1:500-1:5000 IHC: 1:50-1:200

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation, FC = Flow cytometry

### Immunogen information

Immunogen:	Ag4793
GenBank accession number:	BC026335
Gene ID (NCBI):	10391
Full name:	Coronin, actin binding protein, 2B

### Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3. Store at -20°C.