

## STXBP4 Polyclonal Antibody

Catalog number: 14231-1-AP

Size: 39 µg/150 µl


Source: Rabbit

Isotype: IgG


Synonyms:

STXBP4; FLJ16496; MGC149829;

MGC50337; Synip

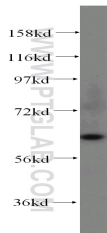

human heart tissue were subjected to SDS PAGE followed by western blot with 14231-1-AP(STXBP4 Antibody) at dilution of 1:1000

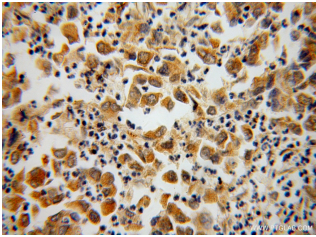

Immunohistochemistry of paraffin-embedded human pancreas cancer using 14231-1-AP(STXBP4 Antibody) at Dilution 1:100 (under 40x lens)

### Applications

Tested applications:	ELISA, WB, IHC
Species specificity:	Human, mouse, rat; other species not tested.
Positive WB detected in:	Human heart tissue, BxPC-3 cells, PC-3 cells
Calculated STXBP4 MW:	62kd
Observed STXBP4 MW:	62kd
Positive IHC detected in:	Human pancreas cancer
Recommended dilution:	WB: 1:500-1:5000 IHC: 1:50-1:200

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation, FC = Flow cytometry

### Immunogen information

Immunogen:	Ag5471
GenBank accession number:	BC041485
Gene ID (NCBI):	252983
Full name:	Syntaxin binding protein 4

### Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3. Store at -20°C.