

NORLAKE®

MINI ROOMS

Controlled
Temperature
Chambers




Mini Rooms are ideal for general-purpose cold storage in laboratories, university and medical research facilities, morgues, agriculture rooms, animal carcass and museum art/collections storage. Mini Rooms are used to store reagents, media, proteins, specimens and other samples and solutions.

SCIENTIFIC

MINI-ROOM WALK-IN CONTROLLED TEMPERATURE CHAMBERS

- Single room & combination cold rooms and freezers
- Temperatures +4°C, -20°C
- Fast and easy installation
- Ceiling mount Capsule Pak™ refrigeration system
- All system components outside allowing more storage inside
- Standard sizes from 4' x 6' to 10' x 14'
- Indoor and outdoor models
- Options: Interior and exterior kick plates (stainless steel or aluminum diamond tread), exterior ramp (for with-floor models), shelving kits, 30" or 36" wide door section upgrade, 24" or 48" LED light fixture (shipped loose)

ALL MODELS SHIP WITHIN
20 BUSINESS DAYS!

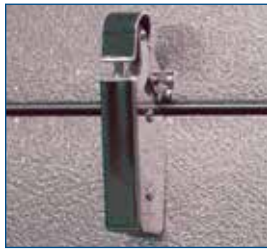

Spring actuated door closer



Adjustable door hinges



LED interior lighting



Deadbolt locking door handle with independent key/padlock feature and inside safety release

CAPSULE PAK™ REFRIGERATION SYSTEMS

- Standard on Mini Room coolers and freezers
- Flush coil makes full use of interior space
- Easy installation with no floor drain required
- All models feature standard cord and plug eliminating the need for field connection
- UL and C-UL electrical listing on complete Capsule Pak refrigeration systems
- AWEF and CARB compliant



Standard LogiTemp™ electronic controller on each Capsule Pak system protects your perishable products while boosting energy efficiency.

NORLAKE®

891 County Road U • Hudson, WI 54016
800-955-5253 • www.norlake.com

RSG
Refrigerated Solutions Group®