

HAAO Polyclonal Antibody

Catalog number: 12791-1-AP


Size: 43 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:

HAAO; 3 HAO, HAAO, HAD,
HAO

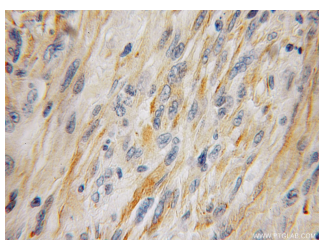

Immunohistochemical of paraffin-embedded human gliomas using 12791-1-AP (HAAO antibody) at dilution of 1:100 (under 10x lens)

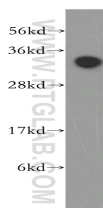

human liver tissue were subjected to SDS PAGE followed by western blot with 12791-1-AP (HAAO antibody) at dilution of 1:300

Applications

Tested applications:	ELISA, WB, IHC
Species specificity:	Human, Mouse, Rat; other species not tested.
Calculated HAAO MW:	286aa, 33 kDa
Observed HAAO MW:	33 kDa
Positive WB detected in	Human liver tissue, mouse liver tissue, rat liver tissue
Positive IHC detected in	Human gliomas tissue
Recommended dilution:	WB: 1:200-1:2000 IHC: 1:20-1:200

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

Immunogen:	Ag3511
GenBank accession number:	BC029510
Gene ID (NCBI):	23498
Full name:	3-hydroxyanthranilate 3,4-dioxygenase

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C.