

NFYA PolyAb Antibody

Catalog number: 12981-1-AP

Size: 74 µg/150 µl


Source: Rabbit

Isotype: IgG


Synonyms:

NFYA; CBF-A; CBF-B; FLJ11236;

HAP2; NF-YA

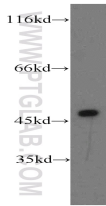

HEK-293 cells were subjected to SDS PAGE followed by western blot with 12981-1-AP(NFYA Antibody) at dilution of 1:500

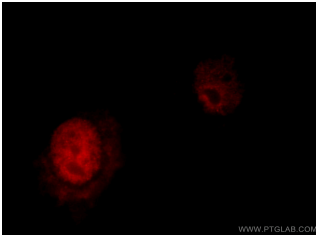

Immunofluorescent analysis of HepG2 cells, using NFYA antibody 12981-1-AP at 1:25 dilution and Rhodamine-labeled goat anti-rabbit IgG (red).

Applications

Tested applications:	ELISA, WB, IF
Species specificity:	Human, Mouse, Rat; other species not tested.
Positive WB detected in:	HEK-293 cells, HeLa cells, HT-1080 cells, Jurkat cells, K-562 cells
Calculated NFYA MW:	318aa, 34kd, 347aa, 39kd
Observed NFYA MW:	40kd, 46kd
Positive IF detected in:	HepG2 cells
Recommended dilution:	WB: 1:200-1:2000 IF: 1:20-1:50

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation, FC = Flow cytometry

Immunogen information

Immunogen:	Ag3604
GenBank accession number:	BC039244
Gene ID (NCBI):	4800
Full name:	Nuclear transcription factor Y, alpha

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3. Store at -20°C.