

## PCYT1B Polyclonal Antibody

Catalog number: 13765-1-AP


Size: 31 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:

PCYT1B; CCT B, CCT beta,  
CCTB, CT B, CTB, PCYT1B

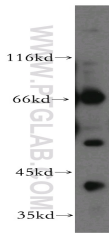

human placenta tissue were subjected to SDS PAGE followed by western blot with 13765-1-AP(PCYT1B antibody) at dilution of 1:500

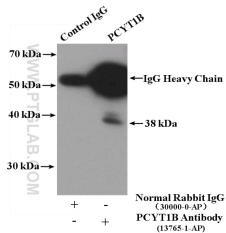

IP Result of anti-PCYT1B (IP:13765-1-AP, 4µg; Detection:13765-1-AP 1:500) with human placenta tissue lysate 4000µg.

### Applications

Tested applications:	ELISA, WB, IP
Species specificity:	Human, Mouse, Rat; other species not tested.
Calculated PCYT1B MW:	351aa, 40 kDa
Observed PCYT1B MW:	42 kDa
Positive WB detected in	Human placenta tissue, human brain tissue
Positive IP detected in	Human placenta tissue
Recommended dilution:	WB: 1:200-1:2000 IP: 1:200-1:2000

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

### Immunogen information

Immunogen:	Ag4748
GenBank accession number:	BC045634
Gene ID (NCBI):	9468
Full name:	Phosphate cytidyltransferase 1, choline, beta

### Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3. Store at -20°C.