

TSEN2 Polyclonal Antibody

Catalog number: 13103-2-AP

Size: 61 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:

TSEN2; HsSen2, PCH2B, SEN2,

SEN2L, tRNA intron

endonuclease Sen2, TSEN2

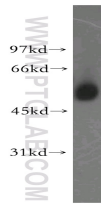

human brain tissue were subjected to SDS PAGE followed by western blot with 13103-2-AP(TSEN2 antibody) at dilution of 1:400

Applications

Tested applications:	ELISA, WB
Cited applications:	WB
Species specificity:	Human,Mouse,Rat; other species not tested.
Cited species:	Human
Calculated TSEN2 MW:	465aa,51 kDa
Observed TSEN2 MW:	51 kDa
Positive WB detected in	Human brain tissue, mouse brain tissue, mouse colon tissue
Recommended dilution:	WB: 1:500-1:5000

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

Immunogen:	Ag3742
GenBank accession number:	BC019582
Gene ID (NCBI):	80746
Full name:	TRNA splicing endonuclease 2 homolog (S. cerevisiae)

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3. Store at -20°C.