

Syntaxin 10 Polyclonal Antibody

Catalog number: 11036-1-AP


Size: 37 µg/150 µl

Source: Rabbit

Isotype: IgG

Synonyms:

STX10; hsyn10, STX10, SYN10,
syntaxin 10

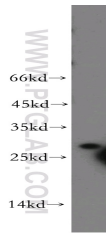

BxPC-3 cells were subjected to SDS PAGE followed by western blot with 11036-1-AP (STX10 antibody) at dilution of 1:200

Background

Syntaxin 10 (STX10), also named as SYN10, belongs to the syntaxin family. Syntaxins are target membrane SNAREs (Soluble NSF attachment protein receptors) which are membrane components that determine the specificity of the docking and fusion of vesicles to target membranes.

Applications

Tested applications:	ELISA, WB
Cited applications:	IF
Species specificity:	Human; other species not tested.
Cited species:	Human
Calculated Syntaxin 10 MW:	28 kDa
Observed Syntaxin 10 MW:	28 kDa
Positive WB detected in	BxPC-3 cells
Recommended dilution:	WB: 1:200-1:1000

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

Immunogen:	Ag1501
GenBank accession number:	BC017237
Gene ID (NCBI):	8677
Full name:	Syntaxin 10

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C.