

TPPP2 Polyclonal Antibody

Catalog number: 13633-1-AP


Size: 27 µg/150 µl

Source: Rabbit

Isotype: IgG

Synonyms:

TPPP2; C14orf8; P18; p25beta

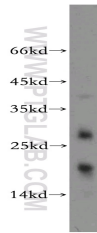

mouse brain tissue were subjected to SDS PAGE followed by western blot with 13633-1-AP(TPPP2 Antibody) at dilution of 1:300

Applications

Tested applications:	ELISA, WB
Species specificity:	Human, mouse, rat; other species not tested.
Positive WB detected in:	Mouse brain tissue
Calculated TPPP2 MW:	170aa,18kd
Observed TPPP2 MW:	18kd
Recommended dilution:	WB: 1:200-1:1000

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation, FC = Flow cytometry

Immunogen information

Immunogen:	Ag4561
GenBank accession number:	BC038970
Gene ID (NCBI):	122664
Full name:	Tubulin polymerization-promoting protein family member 2

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3. Store at -20°C.