

DLEU1 Polyclonal Antibody

Catalog number: 11860-1-AP

Size: 20 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:


DLEU1; BCMS; DLB1; DLEU2;

LEU1; LEU2; MGC22430;

NCRNA00021; XTP6

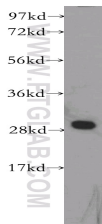

K-562 cells were subjected to SDS PAGE followed by western blot with 11860-1-AP (DLEU1 Antibody) at dilution of 1:600

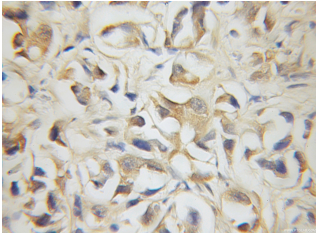

Immunohistochemistry of paraffin-embedded human prostate cancer using 11860-1-AP (DLEU1 Antibody) at Dilution 1:50 (under 40x lens)

Applications

Tested applications:	ELISA, WB, IHC
Species specificity:	Human; other species not tested.
Positive WB detected in:	K-562 cells
Calculated DLEU1 MW:	78aa,9kd
Observed DLEU1 MW:	27-30kd
Positive IHC detected in:	Human prostate cancer
Recommended dilution:	WB: 1:200-1:2000 IHC: 1:20-1:200

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation, FC = Flow cytometry

Immunogen information

Immunogen:	Ag2435
GenBank accession number:	BC020692
Gene ID (NCBI):	10301
Full name:	Deleted in lymphocytic leukemia 1 (non-protein coding)

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C.