

## AGFG2 Polyclonal Antibody

Catalog number: 11919-1-AP


Size: 22 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:


AGFG2; AGFG2, ArfGAP with  
FG repeats 2, HRBL, RAB R,  
RABR

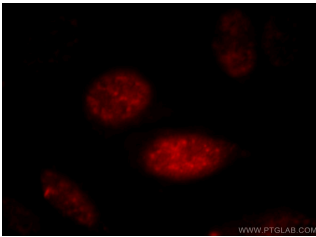

Immunofluorescent analysis of HepG2 cells, using AGFG2 antibody 11919-1-AP at 1:25 dilution and Rhodamine-labeled goat anti-rabbit IgG (red).

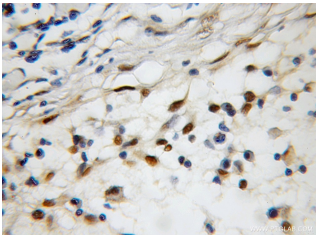

Immunohistochemical of paraffin-embedded human lymphoma using 11919-1-AP(AGFG2 antibody) at dilution of 1:50 (under 10x lens)

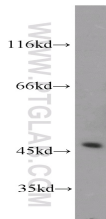

Jurkat cells were subjected to SDS PAGE followed by western blot with 11919-1-AP(AGFG2 antibody) at dilution of 1:800

### Applications

Tested applications:	ELISA, WB, IHC, IF
Species specificity:	Human, Mouse, Rat; other species not tested.
Calculated AGFG2 MW:	481aa, 49 kDa
Observed AGFG2 MW:	45kd
Positive WB detected in:	Jurkat cells, K-562 cells
Positive IHC detected in:	Human lymphoma tissue
Positive IF detected in:	HepG2 cells
Recommended dilution:	WB: 1:200-1:2000 IHC: 1:20-1:200 IF: 1:10-1:100

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

### Immunogen information

Immunogen:	Ag2546
GenBank accession number:	BC017329
Gene ID (NCBI):	3268
Full name:	ArfGAP with FG repeats 2

### Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C.