

CHST1 Polyclonal Antibody

Catalog number: 12504-1-AP


Size: 48 µg/150 µl

Source: Rabbit

Isotype: IgG

Synonyms:

CHST1; C6ST, CHST1, GST 1,
KS6ST, KSGal6ST, KSST

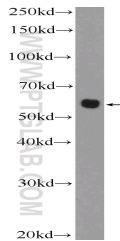

Jurkat cells were subjected to SDS PAGE followed by western blot with 12504-1-AP (CHST1 Antibody) at dilution of 1:600

Applications

Tested applications:	ELISA, WB
Species specificity:	Human, Mouse, Rat; other species not tested.
Calculated CHST1 MW:	411aa, 47 kDa
Observed CHST1 MW:	48-55 kDa
Positive WB detected in	Jurkat cells, mouse brain tissue
Recommended dilution:	WB: 1:200-1:2000

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

Immunogen:	Ag3183
GenBank accession number:	BC022567
Gene ID (NCBI):	8534
Full name:	Carbohydrate (keratan sulfate Gal-6) sulfotransferase 1

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C.