

PCCB Polyclonal Antibody

Catalog number: 11139-1-AP


Size: 45 µg/150 µl

Source: Rabbit

Isotype: IgG

Synonyms:

PCCB; DKFZp451E113, PCCase
subunit beta, PCCB

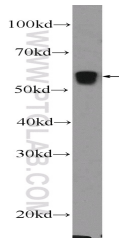

mouse heart tissue were subjected to SDS PAGE followed by western blot with 11139-1-AP (PCCB Antibody) at dilution of 1:600

Applications

| | |
|-------------------------|--|
| Tested applications: | ELISA, WB |
| Species specificity: | Human,Mouse,Rat; other species not tested. |
| Calculated PCCB MW: | 58 kDa |
| Observed PCCB MW: | 58 kDa |
| Positive WB detected in | Mouse heart tissue |
| Recommended dilution: | WB: 1:200-1:2000 |

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

| | |
|---------------------------|--|
| Immunogen: | Ag1608 |
| GenBank accession number: | BC013768 |
| Gene ID (NCBI): | 5096 |
| Full name: | Propionyl Coenzyme A carboxylase, beta polypeptide |

Product information

| | |
|----------------------|---|
| Purification method: | Antigen affinity purification |
| Storage: | PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C. |