

## XPO7 Polyclonal Antibody

Catalog number: 12980-1-AP


Size: 42 µg/150 µl

Source: Rabbit

Isotype: IgG

Synonyms:

XPO7; XPO7, exportin 7, Ran binding protein 16, RANBP16, KIAA0745

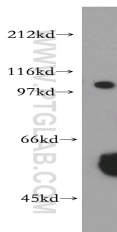

K-562 cells were subjected to SDS PAGE followed by western blot with 12980-1-AP (XPO7 antibody) at dilution of 1:300

### Applications

Tested applications:	ELISA, WB
Species specificity:	Human, mouse, rat; other species not tested.
Calculated XPO7 MW:	1087aa, 124kd
Observed XPO7 MW:	100-110kd
Positive WB detected in:	K-562 cells, human brain tissue, human testis tissue
Recommended dilution:	WB: 1:500 - 1:2000

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

### Immunogen information

Immunogen:	Ag3603
GenBank accession number:	BC030785
Gene ID (NCBI):	23039
Full name:	Exportin 7

### Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3. Store at -20°C.