




- Available with screw cap for low speeds or leak-proof sealing cap for high speeds
- Sealing Cap designed with built-in grips on caps
- Available in autoclavable polycarbonate or polypropylene
- Made in the USA 

Warnings and Guidelines


- May be autoclaved to 121°C at 15 psi for 20 minutes
- Bottles must be filled to 80% capacity for centrifuge work
- Prolonged autoclaving and reuse may cause cracks or discoloration

Foxx Life Sciences' **Autofil®** Centrifuge Bottles are designed specifically to assist Life Science researchers in processing samples for cell suspensions, cell lysates, tissue homogenates, and organ homogenates.

Autofil® Centrifuge Bottles are available in 250ml or 500ml and fit a wide range of rotors. Material options include polycarbonate or polypropylene with your choice of a traditional screw cap or a leak-proof sealing cap with silicone O-ring.

Our centrifuge bottles feature thick walls that hold up to reuse along with wide mouths that simplify filling and pouring liquids. Each **Autofil®** Centrifuge Bottle also features large, clear-molded graduation markings and a frosted area for convenient labeling.

- 1 Screw on or leak-proof sealing cap with built-in grips (leak-proof sealing cap shown)
- 2 Available in autoclavable polycarbonate, chemical-resistant polypropylene, crystal clear polycarbonate, or opaque polypropylene
- 3 High visibility graduations

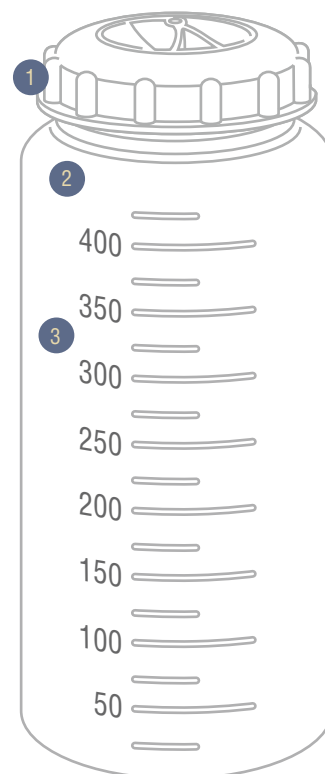



	Item Number	Description	Dimensions	Closure Type	Max. RCF (g) w/cap	Units/ Pack	Packs/ Case
Polypropylene	1801-RLS	250ml PP Bottle with Screw Cap, nonsterile	2.32 x 5.01 in 59 x 129 mm	Screw	8,250	4	9pk
	1803-RLS	250ml PP Bottle with Seal Cap, nonsterile	2.32 x 5.02 in 59 x 129 mm	Seal	27,500	4	9pk
	1805-RLS	500ml PP Bottle with Screw Cap, nonsterile	2.68 x 6.61 in 59 x 129 mm	Screw	8,250	4	6pk
	1807-RLS	500ml PP Bottle with Seal Cap, nonsterile	2.68 x 6.61 in 59 x 129 mm	Seal	13,500	4	6pk
Polycarbonate	1802-RLS	250ml PC Bottle with Screw Cap, nonsterile	2.32 x 5.01 in 59 x 129 mm	Screw	8,250	4	9pk
	1804-RLS	250ml PC Bottle with Seal Cap, nonsterile	2.32 x 5.02 in 59 x 129 mm	Seal	31,700	4	9pk
	1806-RLS	500ml PC Bottle with Screw Cap, nonsterile	2.68 x 6.61 in 59 x 129 mm	Screw	8,250	4	6pk
	1808-RLS	500ml PC Bottle with Seal Cap, nonsterile	2.68 x 6.61 in 59 x 129 mm	Seal	13,850	4	6pk



We recommend that you test bottles for your specific application before use.