

PFDN1 Polyclonal Antibody

Catalog number: 11033-2-AP


Size: 20 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:

PFDN1; PDF; PFD1

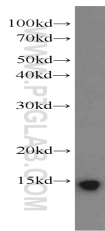

human brain tissue were subjected to SDS PAGE followed by western blot with 11033-2-AP(PFDN1 Antibody) at dilution of 1:600

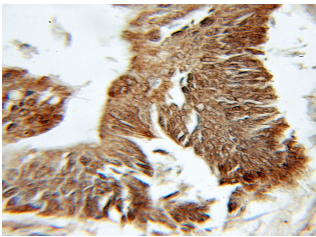

Immunohistochemistry of paraffin-embedded human colon cancer using 11033-2-AP(PFDN1 Antibody) at Dilution 1:50 (under 10x lens)

Applications

Tested applications:	ELISA, WB, IHC
Species specificity:	Human; other species not tested.
Positive WB detected in:	Human brain tissue, human colon tissue
Calculated PFDN1 MW:	14kd
Observed PFDN1 MW:	14 kDa
Positive IHC detected in:	Human colon cancer
Recommended dilution:	WB: 1:200-1:2000 IHC: 1:20-1:200

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation, FC = Flow cytometry

Immunogen information

Immunogen:	Ag1506
GenBank accession number:	BC011869
Gene ID (NCBI):	5201
Full name:	Prefoldin subunit 1

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C.