

SLAIN2 Polyclonal Antibody

Catalog number: 13908-1-AP


Size: 47 µg/150 µl

Source: Rabbit

Isotype: IgG

Synonyms:

SLAIN2; FLJ21611, KIAA1458,
SLAIN motif family, member 2,
SLAIN2

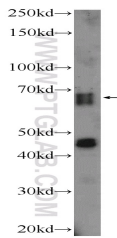

mouse brain tissue were subjected to SDS PAGE followed by western blot with 13908-1-AP(SLAIN2 Antibody) at dilution of 1:300

Applications

Tested applications:	ELISA, WB
Species specificity:	Human,Mouse,Rat; other species not tested.
Calculated SLAIN2 MW:	581aa,63 kDa
Observed SLAIN2 MW:	
Positive WB detected in	Mouse brain tissue, Transfected HEK-293 cells
Recommended dilution:	WB: 1:200-1:1000

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

Immunogen:	Ag4888
GenBank accession number:	BC040993
Gene ID (NCBI):	57606
Full name:	SLAIN motif family, member 2

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C.