

STK31 Polyclonal Antibody

Catalog number: 14140-1-AP


Size: 53 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:

STK31; FLJ16102,
serine/threonine kinase 31,
SGK396, STK31, Sugen kinase
396, TDRD8

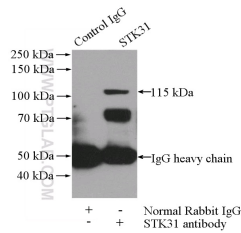

IP Result of anti-STK31 (IP:14140-1-AP, 4µg; Detection:14140-1-AP 1:1000) with mouse testis tissue lysate 2400µg.

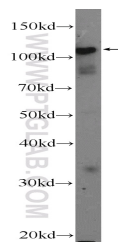

mouse testis tissue were subjected to SDS PAGE followed by western blot with 14140-1-AP(STK31 Antibody) at dilution of 1:600

Applications

Tested applications:	ELISA, WB, IP
Species specificity:	Human,Mouse,Rat; other species not tested.
Calculated STK31 MW:	1019aa,116 kDa
Observed STK31 MW:	115 kDa
Positive WB detected in	Mouse testis tissue, rat testis tissue
Positive IP detected in	Mouse testis tissue
Recommended dilution:	WB: 1:200-1:2000 IP: 1:500-1:5000

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

Immunogen:	Ag5221
GenBank accession number:	BC059374
Gene ID (NCBI):	56164
Full name:	Serine/threonine kinase 31

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C.