

IL3RA Polyclonal Antibody

Catalog number: 13655-1-AP


Size: 20 µg/150 µl

Source: Rabbit

Isotype: IgG

Synonyms:

IL3RA; CD123, hIL 3Ra, IL 3
receptor subunit alpha, IL 3R
alpha, IL 3R subunit alpha, IL
3RA, IL3R, IL3RA, IL3RAY,
IL3RX, IL3RY

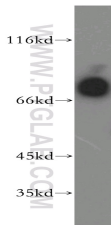

HL-60 cells were subjected to SDS
PAGE followed by western blot with
13655-1-AP(IL3RA antibody) at
dilution of 1:300

Applications

Tested applications:	ELISA, WB
Species specificity:	Human; other species not tested.
Calculated IL3RA MW:	378aa,43 kDa
Observed IL3RA MW:	70 kDa
Positive WB detected in:	HL-60 cells
Recommended dilution:	WB: 1:500-1:5000

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

Immunogen:	Ag4452
GenBank accession number:	BC035407
Gene ID (NCBI):	3563
Full name:	Interleukin 3 receptor, alpha (low affinity)

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3. Store at -20°C.