

POF1B Polyclonal Antibody

Catalog number: 11398-1-AP

Size: 29 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:

POF1B; FLJ22792, POF, POF1B,

premature ovarian failure, 1B,

Protein POF1B

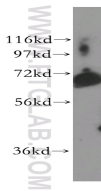

human bladder tissue were subjected to SDS PAGE followed by western blot with 11398-1-AP(POF1B antibody) at dilution of 1:500

Background

POF1B plays a key role in the organization of epithelial monolayers by regulating the actin cytoskeleton. POF1B is a candidate gene for premature ovarian failure (POF); it is mainly expressed in polarised epithelial tissues.

Applications

Tested applications: ELISA, WB

Species specificity: Human, Mouse, Rat; other species not tested.

Calculated POF1B MW: 595aa, 69 kDa

Observed POF1B MW: 69kd

Positive WB detected in Human bladder tissue, mouse bladder tissue, mouse lung tissue

Recommended dilution: WB: 1:500-1:5000

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

Immunogen: Ag1945

GenBank accession number: BC017500

Gene ID (NCBI): 79983

Full name: Premature ovarian failure, 1B

Product information

Purification method: Antigen affinity purification

Storage: PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C.