

RARRES3 Polyclonal Antibody

Catalog number: 12065-1-AP

Size: 56 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:


RARRES3; HRASLS4, PLA1/2 3,

RAR responsive protein TIG3,

RARRES3, RIG1, TIG3

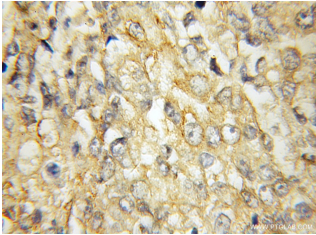

Immunohistochemical of paraffin-embedded human lung cancer using 12065-1-AP (RARRES3 antibody) at dilution of 1:50 (under 10x lens)

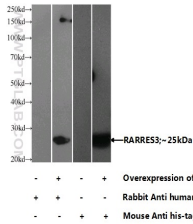

Transfected HEK-293 cells were subjected to SDS PAGE followed by western blot with 12065-1-AP (RARRES3 Antibody) at dilution of 1:1000

Applications

Tested applications:	ELISA, IHC, WB
Species specificity:	Human; other species not tested.
Calculated RARRES3 MW:	164aa, 18 kDa
Observed RARRES3 MW:	
Positive WB detected in	Transfected HEK-293 cells
Positive IHC detected in	Human lung cancer tissue
Recommended dilution:	WB: 1:500-1:5000 IHC: 1:20-1:200

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

Immunogen:	Ag2702
GenBank accession number:	BC009678
Gene ID (NCBI):	5920
Full name:	Retinoic acid receptor responder (tazarotene induced) 3

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C.